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Moc: 0.65W-1.8W

SONIC TOOTHBRUSH FOR KIDS INSTRUCTIONS

IMPORTANT MATTER

Danger

- Do not clean the charging base while power is on.
- Do not immerse the charger or toothbrush for long periods of time in water or other liquids.
- After the charger has finished cleaning, make sure the charger is completely dry and then connects the power supply.
- Check the brush for cracks on a regular basis. Rupture of the brush head can cause a choking hazard.

Warning

- The wires can not be replaceable and discard the charger if the wires are damaged.
- Be sure to replace the original type of charger to avoid danger.
- If the product has any damage to the machine (brush, toothbrush grip and/or charger), please stop using this product.
- Charger can no longer be used near high heat surface.
- Do not let children use this product as a toy.

Caveat


- If your mouth or gums have undergone surgery in the past 2 months, consult your dentist before using this product.
- Consult your dentist if you use this product again with bleeding gums and the bleeding condition lasts for about a week.
- If you have concerns about medical care, consult your doctor before using this product.
- ORO-MED electric toothbrush conforms to the safety standard of electromagnetic equipment. If your body is equipped with a rhythm adjuster or has an implanted device, consult a dentist or an implant device before use.
- This product is only used to clean teeth, gums and tongue, do not use for other purposes. If you feel any discomfort or pain, please stop using it immediately and ask your professional doctor for help.

SAFETY PRECAUTIONS

To reduce the risk of malfunction, electric shock, injury, fire, death, and equipment or property, be sure to observe the following safety precautions.

- Do not disassemble any product other than this product or may cause a fire, injury or electric shock.
- Do not store in infants and young children within the scope of contact. Do not use by infants and young children or may cause an accident or injury due to inadvertent inhalation of moving parts.
- If you do not use this product for a long period of time, unplug the power plug from the home power outlet. Otherwise it may be due to decreased insulation performance caused by leakage, resulting in electric shock or fire.
- People who have severe periodontal disease, have had a history of dental disease, or have a suspected oral condition should consult a doctor before use. Otherwise it may cause tooth loss.
- People with low oral sensitivity should not use this product. Do not force the brush on the teeth or gums, or clean a part for a long time. Do not touch the plastic parts of the unit or the brush head. Otherwise it may cause tooth and gums injury.
- Do not share your brush head with your family or others. Otherwise it may cause infection.
- Do not place metal objects on the charger. Otherwise it will burn because overheating.


PART NAME

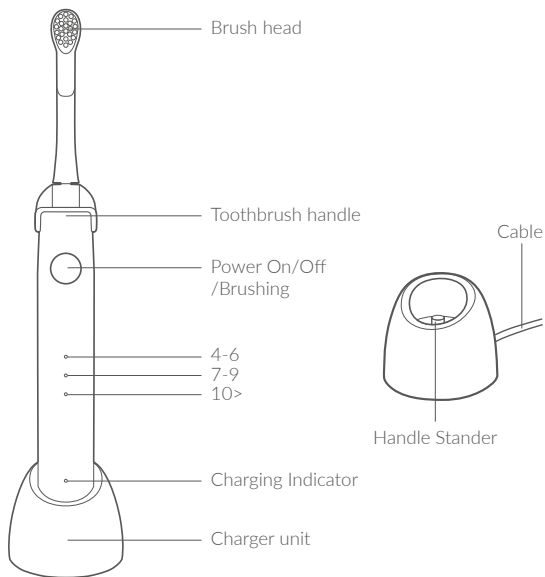
Note: please check the parts and quantity.

INSTRUCTIONS FOR USE


1. Please insert the brush head firmly into the brush head shaft.

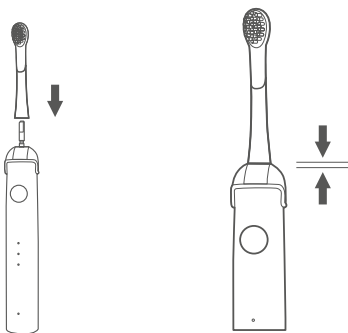
- Install the brush head, first with warm water (do not use too hot water) soaked brush bristles, making the baby brush more comfortable.

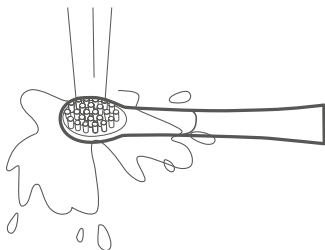2. Please squeeze the toothpaste in appropriate amount.
 - The vertical alignment of toothpaste brush, squeeze proper toothpaste.

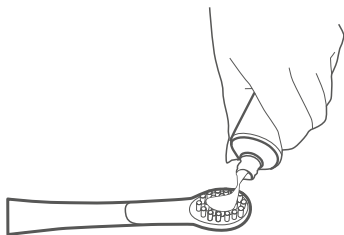3. Press the power button to switch brushing modes.
 - Turn on the power(because its section is dedicated to the baby, micro-vibration, do not worry about the situation of toothpaste splashing), let the baby put it into his mouth and open the automatic brushing mode.

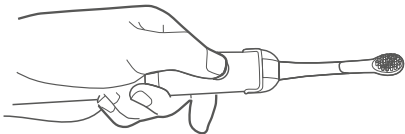4. Turn off the power after brushing.
 - After brushing for 2 minutes, turn off the power automatically and rinse with water.

5. Clean the brush head safely.
- Wash the bristles with clean water, and then brush the head thoroughly and place it in a safe place.

BRUSHING MODES

Professional for children's mouth, 3 brushing modes switch based on children's ages.

- **1 mode** 4-6 years old: soft massage children's gums, gentle cleaning, care baby oral health.
- **2 mode** 7-9 years old: deep cleaning, to help prevent from tooth decay, dental carries and other oral problems.
- **3 mode** above 10 years old: The dental cavity, dental gingiva, pit and fissure, which are difficult to clean, can be thoroughly cleaned.

FUNCTIONAL DESCRIPTION

3 Modes

- **1 mode** 4-6 years old: 28000 strokes/min., gentle swing for two minute.
- **2 mode** 7-9 years old: 28000 strokes/min., gentle clean for two minutes.
- **3 mode** above 10 years old: 28000 strokes/min., soft clean for two minutes.

*This product has intelligent memory function. It will open the last mode automatically when you use next time.

Wireless inductive charging

No need to replace the battery and intelligent wireless inductive charging.

Auto-timer guide

2 min. Intelligent auto-timer, 30 seconds interval timer guide

IPX 7 waterproof


This product can be rinsed by water without any problem.

BRUSH YOUR TEETH TIME

When you use this product, there will be a brief pause every 30 seconds to remind you to change the brushing area.

Clean the various areas of the mouth to achieve a scientific cleaning effect and to develop a good brushing habit.

For your oral health, it is advisable to replace the brush head every 3 months.

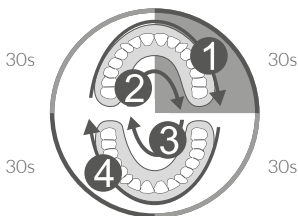
CLEAN THE BRUSH HEAD AND THE HANDLE

Do not let your children clean the toothbrush.


Turn off the toothbrush before rinsing with water.

Advise to change the brush head every 3 months.

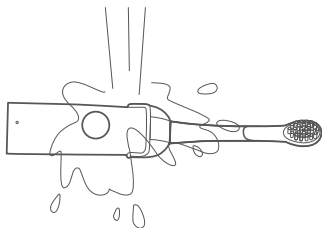
Brush head

- Rinse the brush head with liquid water
- Shake off water after cleaning
- Do not brush the brush head bristles.
- Do not use hot water above 80°C.

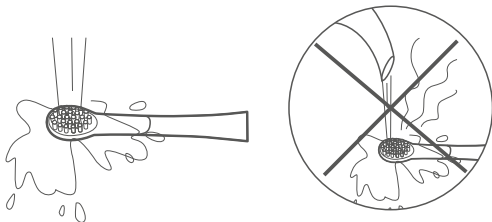Handle

- Rinse the handle with liquid water
- Do not rub or damage the waterproof rubber seal.
- Do not immerse the device into water.
- Do not use hot water above 40°C.
- Wipe the device with dry cloth and put it back to a safe place.

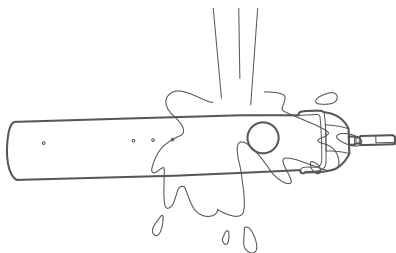Toothbrush charging

- Plug the adaptor into an electrical outlet and put the brush handle onto the charging base, the charging indicator should be lit.
- The full charge needs about 16 hours which enables about 30 days brushing. Power indicator keeps flashing green which means battery low power. Indicator light turns green when charging, and fully charged when indicator keeps flashing green.
- Pull plug out of outlet for saving energy, or leave it and keep the brush on charging base, there is no risk of overcharging of battery.
- If not use for a long time and when you charge it, the power indicator will not be lit immediately which is normal. After pre-charging 30 minutes, the power indicator will be lit.

PARAMETER

Rated input: DC5V 1A
Main working voltage: 3.0V-4.2V
Battery capacity: 600mAh
Power: 0.65W-1.8W

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