PL

EN

DE

FR

ES

PT

IT

CZ

SK

SL

EE

LT

LV

RU

BEAUTIFLY

Shine Ice

PRO
series

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- www.beautifly.eu

- www.youtube.com/beautiflyeu

Light type:	Intense Pulsing Light (IPL)
Adapter:	AC: 100-240V, 50/60Hz +/- 3Hz DC: 12V 6 A 72 W
Product Size:	200 × 170.6 × 87.5 mm
Product weight:	0.5kg
Window size:	4cm2 / 2cm2 / 1cm2
Wavelength:	600nm
Working environment:	5-40°C humidity ≤75%
Storage Environment:	0-60°C humidity ≤75%
Included:	Device, adapter, user manual, safety glasses, razor, cosmetic bag, 2 attachments for different body parts

Product description [page 2- f1]:

1. Place where impulses are released.
2. Ice Cooling area
3. Sensor (skin color and skin contact)
4. Air inlet
5. Air outlet
6. LCD display
7. Mode button
8. Power button
9. Light indicator
10. Adjusting the angle of inclination
11. Pulse emission button
12. Handle
13. Charging port
14. 1 cm2 Filter (for facial hair removal)
15. 2 cm2 Filter (for depilation of armpits and bikini area)

Operation procedure:

To obtain the best results, we suggest that before first use, read the user manual, including all safety rules and warnings, in order to use the device correctly and safely. The B-Shine Ice hair removal device has integrated filters that protect the skin against harmful radiation and burns during the hair removal treatment. The rate of hair growth varies, depending on age, metabolism and other factors. Hair growth is usually divided into three stages:

Growing phase

Hair actively grows along with the bulbs. At this stage, melanin concentration is the highest. Melanin determines the formation of hair color. Melanin concentration is also a decisive factor in the effectiveness of intensive pulsed light technology. Only hair in the growing phase will react to the light hair removal process.

Degeneration phase

Stopping hair growth and rooting the hair bulb.

Reconstruction phase

Mature hair detaches from the bulbs and falls out. The bulb is dormant until the next hair growth cycle.

In the first stage of hair removal with the B-Shine Ice device, light energy is used to act directly on the hair bulbs below the skin. Melanin in the skin and hair bulbs absorb the pulsed light energy emitted by the device. The darker the hair color, the more light energy is absorbed, and this process stimulates the hair follicles to enter the dormant phase. In the second stage, during the depilation period, the hair will fall out naturally, and regeneration will be inhibited, which may take 1-2 weeks. Only hair in the growing phase will be removed. To effectively remove unwanted hair, you should repeat the depilation process once a week for two months.

Display description:

Device Life: When turning on the device or 10 seconds without operation, the device will show the remaining life. Calculation of working times: the light emission times will be shown on the LCD display. After 10 seconds without operation, the device will show the remaining life.

Working mode (6 modes):

1. HR: hair removal,
2. (HR + snowball icon): hair removal + Ice Cooling
3. SR: skin revitalization,
4. (SR + snowball icon): skin revitalization + Ice Cooling
5. AC: acne treatment
6. (AC+ snowball icon): acne treatment + Ice Cooling

Energy regulation:

Total 5 levels, energy intensity changes with different level.

Pulse Mode:

When the display shows one lightning icon, it means manual pulse mode.

When it shows two lightning icons, it means automatic pulse mode.

Countdown on the display after the device has finished working:

Long press the power button to turn off, the screen will show a 30 second countdown. During this period, the fan continues to operate normally to dissipate heat. When the countdown ends, the device is completely turned off.

Control panel specifications:

Power Button:

1. Press the power button for 2 seconds to turn on the device. Press the power button for 2 seconds to turn off the device. The internal cooling system will run for 30 seconds to cool the device after shutdown.
2. Once turned on, the device will automatically select the appropriate energy level according to your skin color when the pulse release point comes into contact with your skin. You can also select the energy level manually. Just tap the power button to adjust the appropriate energy level (level 1-5). It is recommended to start with a low level (level 1), the level can be gradually increased if the skin responds well.

Mode button:

1. Long press **M** button, turn Ice Cooling function on or off
2. Short press **M** button, switch function under HR (hair removal mode), SR (skin revitalization mode) and AC (acne treatment mode).

Broadcast button:

1. Manual Pulse Mode: The display shows the lightning bolt icon, make sure the pulse release point is well attached to the skin (turning on the indicator light). Press the broadcast button once.
2. Auto Pulse Mode: Press and hold the emission button for 2 seconds, the LCD screen shows two lightning icons, make sure that the place of emitting pulses adheres well to the skin. Press the emission button once and the device will automatically start emitting pulses.

Contraindications:

1. The device is not intended for children under 15 years of age. Young people aged 15–18 may use the device with the consent and/or supervision of their parents or guardians.
2. Do not use if the skin in the area to be treated is tanned or red. Using the device on tanned or red skin may cause temporary or even permanent discoloration of the skin.
3. Do not use if you are pregnant, trying to get pregnant, think you may be pregnant or are breastfeeding. The B-Shine Ice device has not been tested for this purpose and therefore should not be used in these cases. Hormonal changes that occur during pregnancy or breastfeeding may increase the risk of skin damage or skin sensitivity.
4. Do not use the device if you suffer from diseases such as epilepsy, epileptic seizures or migraine. Energy flashes can cause severe headaches and migraines.
5. Do not use on areas of skin covered with tattoos, dark brown or black spots (such as dark freckles, birthmarks, or warts). Tattoos or dark spots can absorb too much light energy, which can cause skin swelling, temporary or permanent blisters, discoloration, or skin whitening.
6. Do not use the device if you have a tendency to develop keloid scars.
7. Do not use on the scalp, ears, neck, nipples or mucous membranes (labia minora, vaginal and anal areas). In those parts of the bikini area where the skin is more pigmented, a lower light intensity level should be used. Use in these areas may cause skin swelling, blistering, or temporary discoloration or whitening of the skin.
8. Do not use if the treatment area has any active skin inflammation, including acne, psoriasis, eczema, burns, herpes, abrasions or rash. These factors can make the skin more sensitive to light and make the treatment uncomfortable or cause swelling, blisters or infections on the skin. Once your skin has healed sufficiently, consult your doctor before using the device.
9. Do not use the device if you are taking medications that make your skin sensitive to light, including non-steroidal anti-inflammatory agents (e.g. aspirin, ibuprofen, cetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulfonyleureas, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha Hydroxy Acids (AHAS), Beta-Hydroxi Acids (BHAS), Retin-A, Accutane and/or topical retinoids.
10. Do not use if you have a disease associated with photosensitivity, e.g. porphyria, polymorphic photoeruption, solar urticaria, lupus, etc. These conditions make the skin more sensitive to light and the procedure may be unpleasant or causing skin swelling, blisters or infections.

Handling Warning:

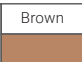
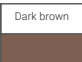
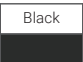
1. Wear safety glasses when operating the device.
2. Clean the area of skin being treated before use. It is not allowed to apply any cosmetics to the skin immediately before the treatment. This may damage the device and affect the effect of the treatment.
3. Adjust your energy level from low to high. The higher the energy level, the better the result. Immediately after the treatment, skin redness may occur.
4. After the treatment, use a moisturizing preparation and preparations with high sun protection.
5. After the treatment, use a damp cotton swab or cloth to clean the lamp from any remaining hair, sebum or sebum.
6. It is recommended that the device and accessories be stored away from water and moisture. Store the device in a ventilated and dry place and keep it away from children.
7. It is prohibited to place metal objects in the heat emission hole. Otherwise there is a risk of electric shock.

Before starting the procedure:

1. Use a hair razor before use.

Note: Unshaven body or uncleaned skin may cause adverse effects of the device:

- Dust in the hair may absorb some of the light energy, which may cause skin pain or redness.
 - The hair on the skin can burn, and in severe cases it will burn the top layer of the skin.
2. Before the treatment, do not use any cosmetics, lotions or deodorants in the place where you are planning the treatment.
 3. Clean the pulse release area with a damp cloth or cotton swab. This is done to remove hair and dust from the treatment window.
 4. If you are using B-Shine Ice for the first time, you can use the automatic pulse mode to decide the energy intensity. You can also Please compare your skin color with the skin color chart to ensure your skin color is suitable for hair removal treatments.

Skin color \ Hair color	Brown	Dark brown	Black	Dark gold	Pale gold	Red	White/gray
							
White	1-5	1-5	1-5	1-5	✗	✗	✗
Beige	1-5	1-5	1-5	1-5	✗	✗	✗
Light brown	1-5	1-5	1-5	1-5	✗	✗	✗
Brown	1-3	1-3	1-3	1-3	✗	✗	✗
Dark brown	✗	✗	✗	✗	✗	✗	✗
Black	✗	✗	✗	✗	✗	✗	✗

5. Start at a lower level and gradually increase until you start to feel pulses on your skin. You should not feel any pain, the procedure should be pleasant.

6. The pulse output area is 4cm². If it is suitable, please use it directly. If you want to perform the treatment on smaller skin areas, such as armpits, the area above the upper lip, or a deep bikini area, you can use a 1cm² and 2cm² overlay.

7. Skin Tone Sensor: Every time the device is turned on, it is locked and cannot be operated. Please place the pulse area close to the skin. When the red sensor light is off, it means the area has passed skin tone detection and it is safe to use the device. After passing the skin tone test, the device will be unlocked and intelligently automatically adjusted to the appropriate energy level (it is recommended to start at a low level when used for the first few uses).

8. If you are using the device for the first time or on a new area of the body, select a small area near the area you intend to treat to perform a skin sensitivity test. Use the lowest device level. Follow the treatment instructions below, but emit no more than 2-3 pulses. After performing the skin test, wait 60 minutes and check if the skin does not react. If after 60 minutes the skin appears normal and there is no redness or other skin reaction, you are ready for the full treatment.

If your device does not detect your skin tone when you turn it on, your device will not be unlocked and you will not be able to operate it.

Product manual:

1. Use a hand razor to shave off unwanted hair.
2. Plug in the power plug, press and hold the power button to turn on the device, press the power button and set the appropriate energy, start from level.
3. Wear glasses to avoid light damage to your eyes. Make sure the light outlet touches your skin, press the emission button. More details can be found in the "Operating Procedure" section.

Series of treatments:

Make sure the pulse outlet touches the skin, press the emission button. The full hair removal cycle lasts 8 weeks. It is normal for hair to continue to grow after the first few uses. A whole series of treatments should be performed.

Treatment frequency:

- (1-4 weeks) - 2 times a week
- (5-8 weeks) - once every 1-2 weeks
- (8 weeks or more) - 1 time in half or one month

1. To effectively remove visible hair and prevent active growth of hair follicles, epilation should be performed at least once every 2 weeks during the first 2 months after the epilation phase.

2. Use it repeatedly according to your hair growth. A person with a hard beard can increase their energy levels and increase their frequency accordingly. Different frequency of use due to different hair growth.

3. Subsequent treatments

* To ensure long-term smoothness of the skin, it is recommended to perform hair removal treatment every 4-6 weeks as a reminder after the series.

Operation procedure:

Hair removal

Application:

- Mustache, armpit hair, body and leg hair, hairline hair on the forehead and bikini area.
- The effect varies depending on skin color. The device is not suitable for white and light colors.

How often to perform hair removal treatment:

Usually it is about 7-9 treatments. Perform the first 3 treatments once a week. For the next 4-9 treatments, perform one therapy every 2-3 weeks. To maintain the effect of hair removal and to prevent the hair follicles from becoming active again, the treatments should be repeated regularly to ensure that the skin remains smooth. If treatments are not continued, the hair follicles will become active again and the hair will begin to grow back. To maintain the results, simply repeat the treatments every 4-8 weeks to ensure that the hair follicles do not become activated and hair regrowth remains inhibited. The intervals between treatments may vary depending on the individual rate of hair regrowth.

How to perform a hair removal treatment:

1. Before using the B-Shine Ice device, prepare your skin by shaving the selected area to be treated. The skin should also be washed.
2. Wear safety glasses to avoid IPL light stimulating your eyes.
3. Select the appropriate intensity (see "select intensity") and start at a low level.
4. Let the device window touch the selected areas vertically and treat the areas one by one. Don't repeat the same place over and over or skip any area.
5. Start on the initial area and repeat step 4, 2-3 times.
6. After the treatment, it is recommended to use a moisturizing preparation after depilation.
7. Turn off the device and use a damp cloth or cotton swab to clean the treatment window.
8. Place the device and accessories in a case or box and store it in a safe place.

Skin photorejuvenation:

Do you want to keep your skin looking young and healthy for years? We use our skin all year round, exposing it to various weather conditions. Exposure to the sun is particularly harmful. Intensive tanning usually leaves marks on the face or neckline, which are often difficult to cover with makeup. Similar discolorations also appear on the backs of the hands over time. How to deal with them?

For people who want to visibly rejuvenate the skin and refresh its appearance, evening out the color and reducing bothersome discolorations, we offer photorejuvenation treatments with the modern B-Shine Ice cosmetic laser. B-Shine Ice laser therapy allows for significant smoothing of the skin and complete removal of age spots and erythema. The aim of rejuvenating selected areas is to use short and safe light pulses directed towards the skin surface. The light wave stimulates skin cells to produce collagen and elastin, which accelerate skin regeneration. The therapeutic effects of IPL radiation include improved blood supply and oxygenation of the skin, an increase in the number of collagen and elastic fibers, better moisture absorption, and regulation of the pigmentation and keratinization processes.

Effects:

- elimination of symptoms of skin aging
- reconstruction of collagen fibers
- firm and elastic skin
- shrinking of pores
- smoothing wrinkles
- improved skin hydration
- lightening discolorations

How often to perform a photorejuvenation treatment:

To obtain the desired effects, a series of treatments must be performed. It is recommended to carry out a series of 8 treatments, one treatment each week.

How to perform a photorejuvenation treatment:

1. Before starting the treatment, you should thoroughly cleanse your skin of makeup and any cosmetics. You should also wipe the lamp.
2. Wear safety glasses to avoid IPL light stimulating your eyes.
3. Select the appropriate intensity (see "select intensity") and start at a low level.
4. Let the device window touch the selected areas vertically and treat the areas one by one. Don't repeat the same place over and over or skip any area.
5. Start on the initial area and repeat step 4, 2-3 times. After the treatment, it is recommended to use a strongly moisturizing preparation.
6. Turn off the device and use a damp cloth or cotton swab to clean the treatment window.
7. Place your device and accessories in a case or box and store them in a safe place.

Acne treatment:

Acne treatment using IPL is intended for people with medium severity of acne lesions. IPL light destroys rapidly multiplying bacteria responsible for inflammatory changes in acne. The action is similar to antibiotics that are used topically or systemically. However, they can cause a number of side effects, damaging the liver and stomach.

B-Shine Ice laser therapy works directly on the causes of acne. It acts on the so-called molecules synthesized by anaerobic bacteria. porphyrins. These make the bacterium particularly sensitive to light. When porphyrin is exposed to visible light, it becomes chemically active and reaches an excited state. This causes the formation of singlet oxygen, which combines with the cell membrane and destroys anaerobic bacteria responsible for the formation of acne lesions and sebaceous gland cells.

IPL light therapy must not be combined with antibiotic therapy. You should also wait a few weeks after completing such treatment.

How often to perform acne treatment:

The treatment period lasts approximately 8 - 16 weeks (2-4 months), and treatments are performed every 1-2 weeks. It is best to perform two irradiations of all affected areas during the procedure, approximately 10 minutes apart. Starting from the second week, the first signs of improvement are visible (the degree of inflammation in the affected areas and a slight reduction in the symptoms themselves occur).

After the entire series of treatments:

Acne spots are almost completely eliminated. However, those that remain - no longer active - disappear over the next weeks.

After healing, it is recommended to use skin photorejuvenation therapy, which is able to stimulate collagen

production in order to reduce the number of acne scars. During the entire treatment, it is recommended to use appropriate care cosmetics intended for oily and acne-prone skin.

How to perform an anti-acne treatment:

1. Before starting the treatment, clean the skin thoroughly from makeup and all cosmetics. You should also wipe the lamp.
2. If inflammation occurs in the area to be treated, the skin should be healed before using the device.
3. Wear safety glasses to avoid IPL light stimulating your eyes.
4. Select the appropriate intensity (see "select intensity") and start at a low level.
5. Let the device window touch the selected areas vertically and treat o bgrays one by one. The interval between pulses should be at least 20 seconds.
6. After the treatment, it is recommended to wash your face and use a moisturizing preparation.
7. Turn off the device and use a damp cloth or cotton swab to clean the treatment window.
8. Place your device and accessories in a case or box and store them in a safe place.

More information about the company and products [page 2 - f2]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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