optimum
med[™]

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Waga: 1,50 kg

NOTICE

- It must be provided in the transformer (DC 12V 2A) or confirm the rated voltage and then use.
- Make sure that products and transformer correctly contact.
- Don't pull the connection by force. The bad connection of transformer will cause the product to stop work.
- Do not use it in the bathroom and humid environments.
- Please do not place the product on the solid surfaces to use, because it will be damaged.
- The product is with heating on the surface, the people who are insensitive must pay attention to it.
- If a product is something wrong, do not privately disassembly it.

WARNING: to reduce fire, burns, electric shock, or other injury risk, please observe the following instructions.

- When nobody cares it, please pull out the plug from the socket.
In addition, when it is out of using or opening it, you must pull out the plug.
- To avoid over heat, electric shock, or injury, do not use the product in Pillows or blankets under operation.
- If the product is close to the children, physical impairment and disability when using, they need comprehensive the guide.
- Please use it according to the manual of the product. Do not use non-manufacturing accessories recommended by other manufacturer.

DANGER:

In order to reduce the risk of electric shock danger, please follow the instructions below:

- After using or before cleaning products, please pull out the plug.
- If the product falls into the water, do not touch the wet product. You must pull out the plug.
- Do not use the product when taking a shower or a bath.
- Do not place the product in water or other liquid.

PARTS NAME AND FUNCTION

ON/OFF

Working time: 15 minutes

It will be closed after using for 15 minutes

DR (DIRECTION)

You can adjust the direction of turning.

It can be used effectively in specific areas. You can massage specific part of your body.

SPEED ADJUSTMENT

You can use it to speed up or slow down, the kneading speed: medium, fast and slow.

HT (HEATER)

Way of using buttons if you want to adjust the buttons. While working, press one time, it will be closed. Press again, it will be open.

REGULAR MASSAGE CAN IMPROVE THE WHOLE BODY FUNCTION.

■ Improve the quality of sleep

Massage can make the body relax, eliminate fatigue, and relieve pressure. It also can improve sleep, and improve sleep quality.

■ Improve body function

Massage can promote blood circulation for the body to provide adequate nutrients, improve lymph detoxification, relax the muscles to maintain body balance of yin and yang.

■ Condition gastrointestinal function

Massage can improve gastrointestinal function, promote the new supersedes the old, improve digestion. Help digestion, relieve constipation and other gastrointestinal adverse circumstances.

MASSAGE OF NECK AND SHOULDER

The new design, exquisite curve design, tightly massage every position.

UNIQUE WAY OF KNEADING MASSAGE

Give up the traditional massage, and have a deep massage on the neck, 4 massage heads knead and massage upward or downward, let you feel relaxed.

THE UNIQUE WARM FUNCTION

Warm function will bring you different feelings. It will bring you the warm feeling while it is having effective massage for you.

HOME USED AND CAR USED

Equipped with external adapter, and you can use it by 12 V vehicle power supply, so you can enjoy massage fun in any place.

USAGE

Be sure the unit is unplugged, cooled and emptied of water. Then store it in a clean and dry place. Before storage just to wrap the whole power cord round the wiring holder. Never hang the unit by the power cord.

HOW TO USE THE MASSAGER

1. Neck massage

2. Leg massage


3. Foot massage

4. Waist massage

5. Abdominal massage

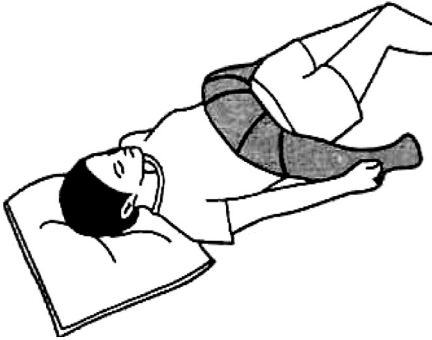CLEANING

- Always unplug the mains adapter from the socket and disconnect it from the device before carrying out any cleaning.
- Clean the device using only the methods specified. Under no circumstances may liquid enter the device or the accessories.
- Small marks can be removed using a cloth or damp sponge, and some liquid delicate laundry detergent where necessary. Do not use any solvent-based cleaning products.
- Do not use the device again until it is completely dry.

MAINS OPERATION:

Input: 100-240 V~; 50/60 Hz

Output: DC 12 V --- 2 A

Rated power: 24 W

ECOLOGY – ENVIRONMENTAL PROTECTION

The “crossed-out trash bin” symbol on electrical equipment or packaging indicates that the device cannot be treated as general household waste and should not be disposed of into containers for such waste.

Obsolete or broken-down electrical device should be delivered to special designated collection points, organized by local public administration, whose purpose is to collect recyclable electrical equipment.

This way, each household helps reduce possible negative effects influencing natural environment and allows the reclaiming of materials which the product is made of.

Weight: 1.5kg 

TROUBLESHOOTING

If you encounter problems when using it, please refer to the following tips and solutions. If the massager is still in the fault state or not available, please turn off the power and unplug the plug.

Users do not disassemble or repair massager by yourself.

You should ask the seller to help you.

Condition	Possible reasons /Solutions
Unable to start the massager	Make sure the power plug is in socket. Confirm the switch is on. Confirm if the button has been pressed.
When using, it suddenly stops to work	<ul style="list-style-type: none"> • Please check whether the power plug is pulled out or not Confirm if the plug has a bad connection. • Maybe press the button control carelessly. • Automatic timer sets for 15 minutes time has passed. • If the overheating protection device starts, please turn off the main power supply, and wait for the massage cooling. If the problem persists, please do not continue to use and contact the company repair service center.