MODEL: MS100K

Please read the instruction leaflet carefully

SPECIFICATIONS

Material: Dupont Nylon/ABS/PETG/TPE Power Supply: 1xAAA battery

Frequency: 16,000 strokes per minute

Product size: 28x25x135 mm

Children should be supervised while using the device. The product is not a toy. Store in place out of children's reach. Children shall not play, clean or repair the appliance. Toothbrush can be used by children and persons with reduced physical, sensory or mental capabilities or lack experience or knowledge if they have been given supervision or instruction concerning the use of the appliance by a person responsible for their safety.

Do not dip the device in water for too long. Do not alter the product in any way other than battery or brushhead change. Do not place the product in high temperatures or expose to direct sunlight. Insert the battery with the correct polarity.

Do not recharge the non-rechargeable batteries. A rechargeable battery must be removed from appliance before being charged. Exhausted battery should be safely disposed of. Do not short-cut any circuits of the appliance.

CHANGING THE REPLACEABLE BRUSH HEAD


To detach: turn brush head anti-clockwise and pull upwise.

To attach: align the new brush head and turn it clockwise.

Note: notch on neck should align with inner nothes of brush head to lock in place.

REPLACING THE BATTERY

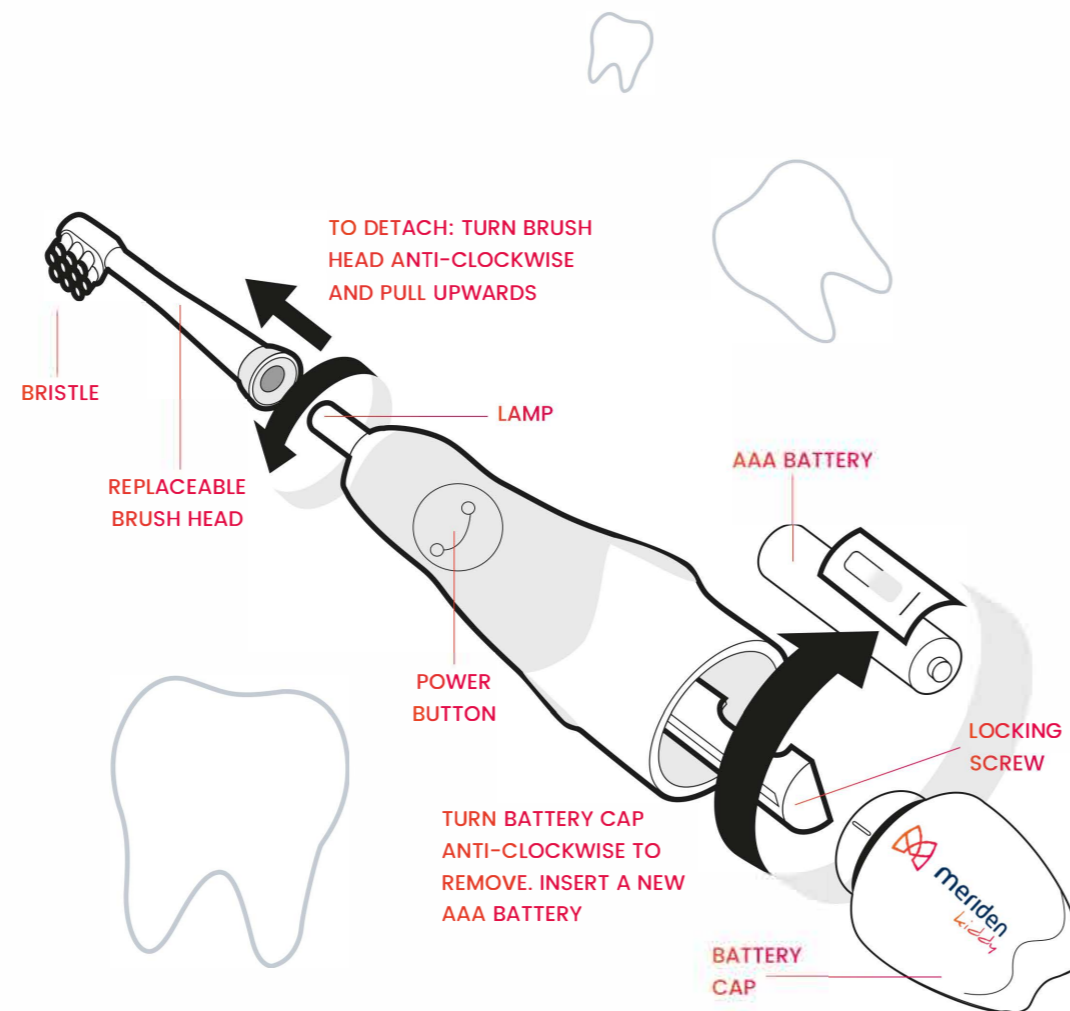- Turn the cap anti-counter to remove brush head
- Using a screwdriver (not included) remove the screw and battery compartment cover
- Insert a new "AAA battery"
- Replace the compartment cover and tighten the screw
- Align line marker on battery cap with the line marker on brush handle. Turn clockwise to replace the battery cap.

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