
SonicBrush Kid G

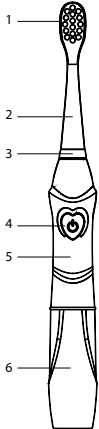EN – USER MANUAL

Box contents

- TrueLife SonicBrush Kid G
- 2x Small brush head + 1x Medium brush head
- Manual and safety instructions

Product description

1. Brush bristles
2. Replaceable brush head
3. LED light
4. On/off button
5. Body of sonic toothbrush
6. Battery cover

Instructions for use


1. Place the brush head onto the body of the toothbrush.
2. Wet the brush head and apply toothpaste to it.
3. Place the toothbrush in the mouth and then switch it on using the on / off button.
4. The LED light will light up, allowing you to view your teeth better.
5. Press the button again to start brushing your teeth.
6. Move the bristles back and forth over your teeth and slowly move between the individual teeth. At the same time, gently push the brush towards and between the teeth.
7. For best results, divide the brushing into 4 sections. Take 30 seconds for each section. After a 30-second cycle, the brush will pause and alert you to move to the next section. After 2 minutes, the brush switches itself off.
8. Rinse the brush head under running water. Wipe the body of the brush with a damp cloth and then wipe it dry with a cloth or towel.

Warning: Children under the age of 7 should not operate the tooth brush without the supervision of a parent or other adult.

Warning: The colour LED indicates that the toothbrush is in use.

Power supply


The toothbrush uses an AAA battery. To insert and replace the battery, proceed as described below.

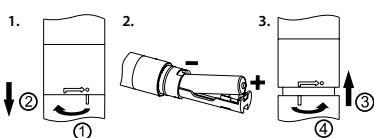1. Remove the cover as shown in figure 1.
2. Insert a new AAA battery or replace the AAA battery as shown in Figure 2.
3. Replace the cover as shown in Figure 3.

Replacing the brush head

To replace the brush head, proceed as described below.

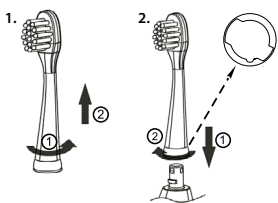1. Rotate the brush head and pull it upwards as shown in figure 1.
2. Insert the new brush head as shown in figure 2.

Warning: The notch on the inside of the brush head should align with the internal notches on the toothbrush body so that the brush head snaps into place.

Maintenance

Toothbrush:

1. Remove the brush head.
2. Wipe any dirt off the surfaces with a soft cloth. Do not immerse the whole body of the brush in water.
3. Rinse the brush head and its bristles after each use.
4. Do not wash the toothbrush or any of its components in the dishwasher.

Specifications

Material	ABS/TPE/ DuPont Nylon
Power supply	1x AAA battery
Battery life	150 days (alkaline batteries, used twice a day for 2 minutes)
Oscillation frequency	24 000/min
Dimensions	19x19x188 mm
Weight	39 g (toothbrush with battery and brush head)

Product liability:

TrueLife products are guaranteed to be free from manufacturing defects. In the event of such defects, elem6 s.r.o. will repair or replace the defective part or product with a comparable part or product at its sole discretion. The warranty will be deemed void if the device is altered or tampered with. The warranty does not cover damage resulting from use, accident or normal wear and tear. elem6 s.r.o. assumes no responsibility for any accident, injury, death, loss, or other claims related to or resulting from the use of this product. elem6 s.r.o. shall in no event be liable for incidental or consequential damages in connection with or arising from the use of this product or any of its components. Return or exchange of spare parts and / or products may be subject to shipping, handling, and replacement or replenishment fees.

Printing errors and changes of information reserved.

SAFETY GUIDELINES WHEN OPERATING THE DEVICE:

1. Read the leaflet included in the packaging carefully.
2. The latest version of the manual, including instructions for operating the device, can be downloaded at: www.truelife.eu.
3. Printing errors and changes of information reserved.

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