TAKE CARE!

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GB	Electronic scale	2
PL	Waga elektroniczna	3
BG	Електронна везна	4
CZ	Digitální váha	5
DE	Elektronische Wiege	6
ES	Báscula electrónica	7
FR	Balance électronique	8
GE	ელექტრონული სასწორი	9
GR	Ηλεκτρονική ζυγαριά	10
HR	Elektronska vaga	11
HU	Digitális babamérleg	12
LT	Elektroninės svarstyklės	13
LV	Elektroniskie svāri	14
NL	Elektronische weegschaal	15
PT	Balança eletrónica	16
RO	Cântar electric	17
RU	Весы электронные	18
SE	Elektronisk våg	19
SK	Digitálna váha	20
TR	Elektronik tartı	21
UA	Весы электронные	22
XS	Elektronska vaga	23

Read the user's manual and warnings before using the device for the first time. Pay special attention to safety tips to avoid accidents and/or damage to the product when using it. Keep the user's manual for reference purposes!

The modern electronic scales are a device for checking the baby's weight. The device is especially useful for weighing infants and premature babies, where the weight is an exceptionally important development factor.

With your needs in mind, we created a product that is easy to use, accurate and highly reliable. The scales have a taring function across the entire weighing range. The scales automatically lock the result for easier reading. We used our best efforts to make sure the product is of top quality, comfortable and safe during weighing.

The weight limit (max. 20 kg) is designed for babies.

TECHNICAL DATA


- Maximal load of the tray - 20 kg
- Minimal load - 0.1 kg
- Weighing accuracy - 0.01 kg
- Measuring unit - kg / lb
- Automatic deactivation after 25 seconds
- Hold function (easier weighing of moving babies - the result is automatically locked)
- Taring function
- Power - 2 x 1.5 V AAA batteries
- Large 60 x 21 mm LCD display
- Dimensions - 550 x 330 x 27 mm

SAFETY PRECAUTIONS

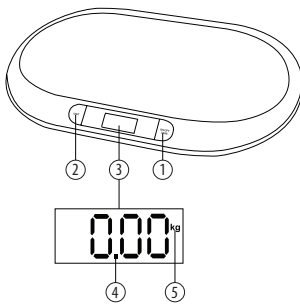- Use to the product for its intended purpose.
- Do not overload the scales.
- The maximal load is 20 kg.
- Avoid shaking or hitting the scales.
- Do not submerge in water.
- When in storage, do not place anything on the scales.
- Do not use or store the scales in humid environment.
- The optimal temperature for operating and storing the scales is from 5°C to 30°C.
- When in use, the scales should be placed on flat, dry and hard surface.

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