# truelife

## RelaxFeet F3


### EN – INSTRUCTIONS

#### PACKAGE CONTENTS

- Massage device RelaxFeet F3
- Power adapter with cable
- Instructions
- Safety instructions

#### PRODUCT DESCRIPTION

1. Removable and washable covers
2. Two heating rotary heads
3. Non-slip mat
4. On / off button
5. Button to change direction of rotation / turn off rotation
6. Heating button
7. Product power cable
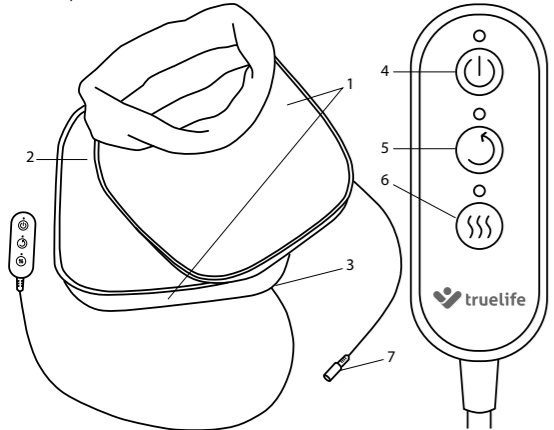
#### USE

Connect the power plug to the product's power cord and the adapter to the outlet.

Press the button (1) to turn on the device. By pressing the button for the first time (2) change the direction of rotation and press it again to turn off rotation. By pressing the button for the first time (3) turn the heater off and press it again to turn the heater on.


You can use the product either for a foot massage, or after removing the top cover for an upper or lower back massage. The massager has a built-in timer and turns off automatically after 15 minutes.

#### Note:

- Use the device no more than 4 times a day. Do not use the device for more than 30 minutes at a time. Wait at least 30 minutes between next use.
- In case of any health complications, consult your doctor about the use of the product.
- Do not use the product on places where you feel pain or have injured skin.
- Never lie down, sit or stand on the massager with your full weight. Use the massager only while sitting.

#### Foot massage

1. Leave or remove the top cover on the massager. By leaving the top part of the sleeve, more heat is retained in the pocket when using the heater.
2. Sit on a chair or sofa and place your feet on the upper part of the massager.
3. Turn on the massager and set the direction of rotation and heat according to your preferences.
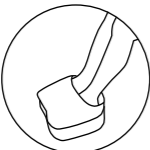
#### Massage of the back

1. Detach the upper part of the cover of the massager from the lower part.
2. Sit on a chair or sofa and place the massager between your back and the backrest.
3. Turn on the massager and set the direction of rotation and heat according to your preferences.

#### MAINTENANCE

Unplug the device from the mains before cleaning. The upper and lower cover are machine washable. Unzip the sleeves and machine wash at 30°C, gentle wash (do not bleach, tumble dry or iron). Allow the sleeves to dry before returning the sleeves back to the device.

#### SPECIFICATIONS

|                         |                                      |
|-------------------------|--------------------------------------|
| Heating                 | Yes (40 °C)                          |
| Spinning back and forth | Yes                                  |
| Timer                   | 15 minutes                           |
| Number of massage heads | 2 (each head has four massage nodes) |
| Controller              | Yes                                  |
| Power cord length       | 3 m                                  |
| Input voltage           | 100–240 V / 50Hz/60Hz                |
| Output voltage          | 12 V – 2 A, 24 W                     |
| Dimensions              | 300 × 300 × 200 mm                   |
| Mass                    | 1209 g                               |

#### MANUFACTURER'S LIABILITY

TrueLife products are warranted against manufacturing defects. Obligation of elem6 s.r.o. in the case of such defects, it is under warranty to repair or replace the defective part or product with a comparable part or product at its own discretion. Remodeling, alteration or other unauthorized handling of the product means exceeding the scope of this warranty. This is even if the product is lost or damaged as a result of a defect, negligence or other error. The warranty does not cover damage resulting from use, accident or normal wear and tear. The company elem6 s.r.o. assumes no responsibility for any accident, injury, death, loss or other claim related to the product or resulting from the use of this product. The company elem6 s.r.o. shall in no event be liable for incidental or consequential damages related to or arising out of the use of this product or any of its components. Returns or exchanges of spare parts and/or products may be subject to shipping, handling, exchange and/or restocking fees.

#### SAFETY WHEN WORKING WITH THE DEVICE

1. Read the package safety information carefully.
2. For proper handling of the product, read the current version of the manual, which is available for download on the page: [www.truelife.eu](http://www.truelife.eu).

Typographical errors and information changes reserved.

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