
BEAUTIFLY

B-Bubble

Body
series

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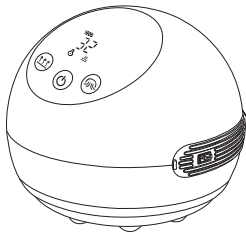
Device description

- . Temperature mode indicator
- . Temperature regulation
- . Power button - on / off
- . Vacuum mode indicator
- . Vacuum adjustment
- . Non-slip silicone
- . Charging interface
- . 8 massage balls

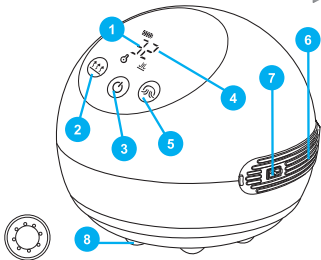Technical data B-Bubble

Size	80*75mm
Net weight	235g
Rated voltage	3.7V
Working current	2A
Power	4W
Red light temperature	36-45°C
Number of gears	5

Red light

It is used to reduce laxity and thicken all skin types. The wavelength is $625 \text{ nm} \pm 3 \text{ nm}$, which allows deep drainage of all layers of the skin. Bio-effects: stimulating cells to produce collagen and elastin faster, increasing cell viability, accelerating blood circulation. It has a strong firming effect, thanks to which it perfectly tightens and tones the skin. It leaves a lifting effect and a smooth body.

Vacuum massage

It uses negative pressure, thanks to which it sucks the skin and goes to the connective tissue. It allows non-invasive cosmetic and rehabilitation effects in deep skin and muscle structures. It also affects the circulation and fat cells. Bio-effects: stimulates the muscles and skin to regenerate, improves circulation, helps to reduce the amount of adipose tissue, has a relaxing and detoxifying effect.


Magnetotherapy

The energy magnet regulates the magnetic balance in the body. It cleanses the blood, improves its flow throughout the body and vigorously stimulates internal organs.

Device operation

Before using the device, use body oil. The care oil acts directly on the skin and helps to enhance the effects of treatments. Use the right amount of oil, depending on the size of the treatment area. Before the procedure, place the oil on your hand, then apply it on the skin and massage it in to absorb.

Attention! Before using the machine, insert the cotton filter into the filter cap, and then insert the filter cap into the suction opening of the device. Thanks to this, you will effectively block the entry of oils into the device, which will extend its life (after several uses, you can wash the filter cotton, then gently wipe it with a paper towel and let it dry, which protects the device better).

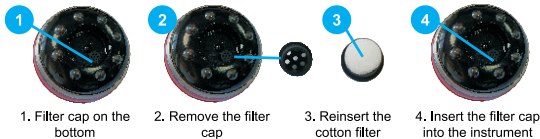Use of functionality:

01 / A long press of the button [figure 1] starts the device which goes into standby mode and the screen shows 00.

02 / Long press of the button [figure 1] turns off the device.

03 / Briefly press the button [figure 2], the vacuum icon [figure 3] will appear, the machine will go into the vacuum 1 mode. (-0 shutdown -1 recycling), when the gear is at 0, the vacuum is turned off.

04 / Pressing the button shortly [figure 4], the temperature icon [figure 5] appears, there is a red light at the bottom, a short press of the button [figure 4] allows you to change gear (0-1-2-3-4-5-0 shutdown-1 recycling). When the device is in gear 0, the heating is turned off with red light.

05 /

- a. In standby mode, the machine will turn off automatically after 5 minutes of inactivity.
- b. Vacuum mode, intensity adjustment: 1st gear - 2nd gear - 3rd gear - 4th gear - 5th gear
- c. The temperature mode corresponds to the red light, and the intensity control is: 1st gear - 2nd gear - 3rd gear - 4th gear - 5th gear
 temperature, 1st gear $36^{\circ}\text{C} \pm 2$
 temperature of the 2nd gear $38^{\circ}\text{C} \pm 2$
 3rd gear temperature $40^{\circ}\text{C} \pm 2$
 temperature of the 4th gear $42^{\circ}\text{C} \pm 2$
 5th gear temperature $45^{\circ}\text{C} \pm 2$

06 / The device stops working when it is connected to the USB port to start charging.
 The device cannot work while charging.

Note 1: Do not use the device when the body is wet.

Note 2: the time of a single treatment should not be too long. To avoid discomfort, it is recommended that the treatment lasts about 5 minutes.

Ways of performing treatments

During treatments, movements in the appropriate direction, appropriate for each part of the body, should be used. After applying the conditioning oil, use the device according to the following operating methods:

- . Back: Massage your back from top to bottom
- . Feet: Massage from heel to toe
- . Belly / buttocks: massage from top to bottom or make circular movements
- . Upper limbs: move from top to bottom
- . Legs: use a movement from the calf up

Tip: For faster cellulite reduction, anti-cellulite oil is recommended.

Issues that require attention after treatments:

1. The treatments open sweat pores and release pollutants. This will consume some of your body fluids. Therefore, after each treatment, you should drink two glasses of water to replenish the fluids used up in the body, support the metabolism and accelerate the excretion of metabolites. After the procedure, you should rest for a while.
2. Due to the negative pressure, the device may cause skin hypersensitivity. If you feel discomfort, suspend use for a few days, and then use the device again after the treatment symptoms have subsided.
3. The recommended treatment duration is 30-60 seconds at a time.

Note: The use of electronic massage devices may harm people with certain physical illnesses. Therefore, if in doubt, consult your doctor before using the device.

Charging the device

1. The device comes with a lithium battery. When the power of the device is low, the battery symbol on the screen starts flashing. To charge the device, connect the appropriate power cord.
2. Use the DC3.7V power adapter, connect the USB and charging cable provided with the device to the power adapter, and insert the other end into the device's Type_Cjack socket to start charging, the battery symbol on the screen will flash. When the battery symbol stops flashing, charging is complete.
3. Use the DC3.7V, 2A or PC USB power adapter to charge the product.
4. The device is idle while it is being charged and should not be used.

Usage Notes

1. If the charging cable or plug is damaged, please do not disassemble it without authorization, please contact the seller timely for repair and maintenance.
2. Disconnect the plug from the power supply after use or before cleaning the device to avoid damage to the device or personal injury.
3. Do not use non-recommended accessories.
4. Do not use the device with a damaged housing or charging cable.
5. Do not use this device when drunk or feeling unwell.
6. Do not use excessive force when using this device to avoid injury.
7. Do not wash the device with water.
8. Do not use this device in places with high temperature, high humidity, inflammable properties and exposed to electromagnetic radiation.
9. Do not use this device while walking or driving a vehicle.
10. Do not use this device in environments with heavy dust or corrosive gases.
11. The battery used in the device is a dedicated rechargeable battery and cannot be used in other products.
12. Protect the device against mechanical damage and in no case disassemble the device.
13. Do not connect the positive and negative terminals of the battery directly to metal parts.
14. Do not store the product with necklaces, hairpins, etc.
15. Do not charge, use or place the product in places with high temperature, such as a stove or hot sun. The above conditions may cause overheating, fire, cracking, etc.

Note: The device is a personal care product, it cannot replace drugs and medical equipment.

Contraindications:

- varicose veins
- menstruation
- damaged and infected skin surface
- breaking blood vessels

Note: The device must not be used for treatments on the face, neck and chest. In the event of hypoglycaemia and hypertension, the device should be used after consulting a physician. The device should not be used immediately after eating due to the higher blood concentration in the digestive system. If you experience any discomfort or tingling sensation while using the device, please stop using it immediately.

Cleaning and maintenance

1. When cleaning and maintaining the product, turn off the appliance and clean the product with a damp towel. If the product is particularly dirty, it can be disinfected with medical alcohol (75% concentration).
2. During daily maintenance, the device should be wiped with a soft cloth with neutral detergent, do not use aggressive cleaning agents, do not soak or rinse in water.
3. During the cleaning process, prevent water or other liquids from getting inside the device to avoid the risk of short circuit and electric shock.

Storage

The device should be stored out of the reach of children, away from moisture, high temperature, direct sunlight and water splashes. When the device is not used for a long time, turn off the power and recharge it every two months.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out basket on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may lead to health detriment or environmental pollution during further disposal or recycling of the equipment.

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The importer hereby declares that the Beautifly B-Bubble device complies with Directive 2014/35 / EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use in certain voltage limits and Directive 2014/30 / EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility. The full text of the declaration of conformity is available at the following internet address: ([link below](#))

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