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DO ZĘBÓW

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IMPORTANT NOTES

Declaration:

This product has reached the highest relative to the manual toothbrush ten times the plaque removal effect, there may be itching or slight bleeding for the first time using, this is because the gum is not used to efficient clean, this phenomenon will stop in 1-2 weeks, if more than two weeks, it may indicate that your gum abnormality in this case should stop using the toothbrush and consult a dentist.

Warning:

1. Orthodontics, dental repair, dental implants and other users, please use them under the guidance of a dentist.
2. This product is only designed for cleaning teeth, gums and tongue. Do not use it for any other purpose. Stop using this product and contact your doctor if you feel any discomfort or pain.
3. Please replace the toothbrush every three months. If there is any wear and tear, it should be replaced in advance.
4. If you are using the toothpaste containing peroxide. Sodium bicarbonate or other bicarbonate, please clean brush with soap and water after using which can prevent possible plastic cracking.



Prohibitions:

1. Do not immerse the charging seat in water or any liquid.
2. This product contains built-in rechargeable batteries, do not charge it at high temperature.
3. Do not knock or drop this product. please stop using the product due to collision or drop.
4. Do not reinvent or repair this product by yourself.
5. please stay away from strong electromagnetic field (such as induction cooker)in charging.

Notice:

1. Do not share the brush head with somebody else.
2. Do not store the power cord around the adapter in storage.
3. Do not place metal objects in the charger, such as coins or paper clips.
4. If you do not use it for a long time, please unplug it so as not to cause fire or danger.

PRODUCT DESCRIPTION

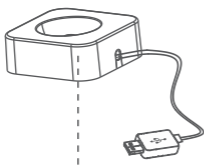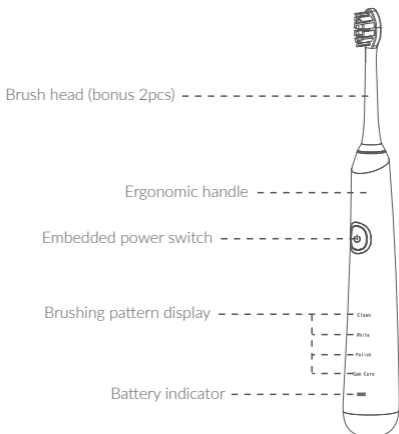
Charger base

User manual


BRUSHING/STEP ONE

Fix the brush head on the toothbrush and make sure there is a certain gap between the two parts. (Gap is to ensure that the toothbrush head does not collide with the main motor during the high frequency vibration, so that the motor has more power).

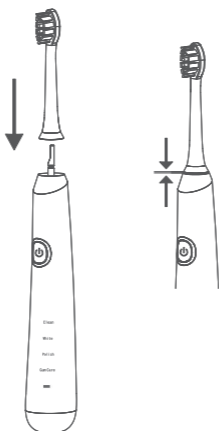
BRUSHING/STEP TWO

Wet the brush head and squeeze right amount toothpaste on brush head.

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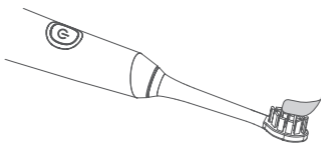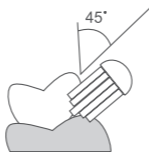Switch on: short press

- **Mode switch:** short press between switch-on condition.

- **Turn-off:** stop working after mode switch between switch-on condition.

Clean

Effective removal of plaque (standard mode).

White

2 minutes whitening pattern to remove surface stains.

Polish

A combination of cleansing and whitening patterns, gentle stimulation and massage of gums, improves gum health.

Gum Care


Extremely soft patterns for sensitive teeth and gums.

CHARGING

Connect the charger base then clean the toothbrush handle before put on the charger base for charging (direction as shown in the figure)

- In the charging process charging indicator light will turn on, the indicator light flashes intermittently

- The indicator light is always on after full charge(full of electricity above sixteen hours);the indicator light flashes quickly when electricity is not enough.

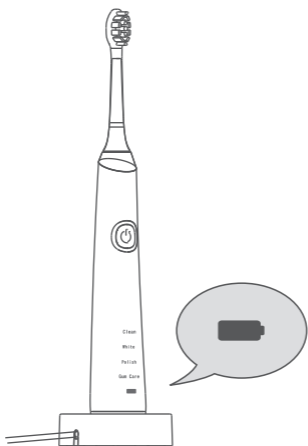The indicator light flashes intermittently in the charging process

Product parameters

Product model	ORO-SONIC PROFESSIONAL
Working pattern	cleaning/whiten/polish/gum care
Working voltage	3.7V
Rated power	3W
Switch way	Tact Switch

Display	LED light
Life time	360minutes
Charging time	16-20hrs
Charging current	50mA
Battery parameter	li-ion battery 2000mAh
Electricity reminding	low electricity indicator
Waterproof class	IPX7
Standard	GB4706.1-2005 GB4705.59-2008

DAILY MAINTENANCE

Brush head

Clean wash by water after each time using ,dry naturally;

Pls use original induction replaceable brush head for your better experience and extended the product life.

Toothbrush handle and charger base

Directly clean by water ,dry naturally;

Pls don't clean with above 40 degree water or acid cleaning liquid.

Induction wireless charging

The toothbrush will be powered up before leaving the factory, in order to maximize the life and efficiency of the built-in lithium battery.

Pls charging for continuous 24 hours for first time using.

Charging indicator light flash continuous shown that low battery power ,pls charging soon.

Please recharge the battery when it is used up, and extend the service life of the lithium battery; The indicator lights flash in the form of breathing in the charging process, and the indicator light is always on after full charge;

When charging, do not put metal objects on the charging base to avoid danger;

At least once a year, the electric toothbrush will be discharged completely.

Do not fall off

Do not use if the body or brush is damaged.

Keep away from heat

The body and charging base shall not be used in a high heat environment.

Please do not repair by yourself

If you need any repairs, please contact our after sales department.

COMMON QUESTIONS

Q: The toothbrush seems to vibrate more and more weakly?

A: 1.Might low battery power ,pls charging soon.
2.Might the brush head connect too tight with main body, pls replace the brush head and fix again.

Q: The battery life time is still short after full charging?

A: 1.pls confirm full charge no less than 20hours each time.
2. pls make sure the toothbrush place vertical-ly.skew place ment will influence the charging effect.

3. if problem still unsolved pls contact the after sale service department for replacement while the battery is out of life.

Q: The charging indicator is not working when the toothbrush was not using for a long time?

A: Pls do not forget to use the toothbrush regularly! If the tooth brush has not been used for more than three months ,the battery will enter into auto-protection pattern.pls continuous full charging above 24hours and check if it could work again. If not pls contact the after sale service department for repair.

Q: What should i do if the toothbrush can not switch on?

A: Pls check if the toothbrush has low battery power and full charging ,if still can not work pls contact after sale service department for repair or replacement.

Q: Could mains cord replaceable??

A: The mains cord can not be replaceable,if the cord broke pls drop it correctly,

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