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MODEL: ORO-SMILE
(MINT/PINK/NAVY BLUE)

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INTRODUCTION

Thank you for choosing OROMED electric toothbrush. Please check if the package is complete before using.

PACKAGE CONTENT

Handle

Toothbrush head

User manual


Type - C cable

ATTACHING THE TOOTHBRUSH HEAD

Set brush head into brush handle and keep same side. Firmly press the heads down on the handle before it stops.

DESCRIPTION

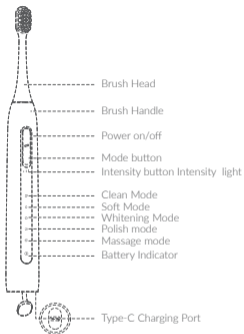MODE SELECTION

Press Power on/off button to power on the toothbrush, press Power on/off again to change mode, the mode will be selected if no further operation in 5 seconds.

INTENSITY SETTING

Power on or Power off the toothbrush, press intensity button to choose 3 different intensity levels: low, medium and high. Intensity levels per mode.

-Low:

1 LED indicator light.

-Medium:

2 LED indicator light.

-High:

3 LED indicator light.

POWER OFF

1. Press Power on/off button about 1.5 seconds while brushing;
2. After all modes are switched in sequence once, press On/Off button once to power off;
3. Automatically power off while charging

CHARGING

1. If the indicator flashes after toothbrush power on, it should be re-charged the indicator will be solid on when the toothbrush fully charged.
2. Wipe the bottom and keep it dry before charging
3. Please keep the bottom cover close after charging

FAQ

1. What should I do if the toothbrush I just received cannot be powered on? It may be caused by insufficient battery power. Please charge it before using.
2. Why the working time of toothbrush less than described? The battery life of more than 90 days is tested by using twice a day and once for two minutes. The working time of toothbrush is subject to the using habits. The battery life is for reference only
3. Why does the vibration decrease after a period of use? As the use time increases, the power decreases, which will cause a decrease in the vibration. You just need to charge it.

IMPORTANT

1. Before discarding the toothbrush, the battery must be removed:
2. The battery must be removed in a power-off condition,
3. The battery should be disposed of safely;

4. Disassembly must be done by professionals only

Steps of battery removal:

1. Remove the bottom cover,
2. Remove the movement support assembly from the body shell;
3. Heat the battery pins with electric ferrochrome to remove the battery.

MAINTENANCE

1. Wash the brush head and bristles after each use and allow them to dry naturally.
2. Do not pull the bristles, otherwise, they will be easily deformed, which may lead to falling off.
3. Do not immerse the handle of the toothbrush in water.
4. Do not clean the toothbrush using water higher than 40°C or acidic cleaning solution.
5. Please clean the bottom cover of the toothbrush with a soft cloth.
6. Wipe the plug with a dry cloth to avoid short circuits,

7. It is recommended to replace the toothbrush head every three months.
8. When you use it for the first time, make sure the battery is exhausted before charging.
9. If you have not used it for more than 6 months, charge it for 3-5 hours before using.
10. During use, if the charging indicator flashes continuously for a few seconds and then stops, it indicates that the battery is low.
11. The indicator flashes slowly when it is being charged, and stays on after it is fully charged. If it has not been charged for a long time, the indicator will light on in a recharging delay of 6-10 minutes, which is normal.

AFTER-SALES POLICY

Thank you for choosing our products. We provide a 2-years free warranty for the product, and unopened toothbrush head .

In case of any fault within 24 months from the date of purchase, if it is verified by OROMED After-sales Service Center as a quality problem, you will be provided with a free repair.

Any damage or defect which is not attributable to the product quality is not covered by the warranty, including but not limited to:

1. Failure or damage caused by acts of God, earthquakes lightning strikes, abnormal voltages, and environmental factors;
2. Dismantling, modifying or replace with non-original parts.
3. Do not operate according to the manual or cause damage due to improper use, such as using the wrong current and voltage, plugging the product into an unqualified outlet or other man-made damage.
4. If you need repair within the scope of the warranty please complete the warranty card. If we are unable to perform the obligation of free warranty service due to inaccurate information, you shall bear the responsibility.
5. For details of service and warranty, please visit official website or send email to serwis@oromed.pl

WARRANTY CARD

Thank you for choosing our product. OROMED will provide you with 24 months of after-sales service. In case of any non-human-made damages during this period, OROMED will provide replacement instead of repair. Please complete the following form and send the faulty product back to the seller.

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