optimum

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INTRODUCTION:

Proper use produces a measurable increase in blood circulation in the skin of the feet which manifests itself in a pleasant and lasting feeling of warmth (rise in skin temperature). Foot massager baths at regular intervals can produce improved blood circulation in the feet and indirectly in the entire leg area as well.

Temperature-controlled full-length heaters under each of foot rests maintain the desired temperature throughout your massager when you set the rotary switch of the unit.

There are about 400 stimu-nodes act as gentle fingers to make contact with the soles of your feet and give you most effective and invigorating massage possible. You can easily control their effectiveness by simply pressing down lightly or letting up the pressure with your feet. There are 70 air jets provide a relaxation warm bubble bath for your feet.

READ INSTRUCTIONS BEFORE USING:

When using electrical appliances, especially when children are present, basic safety precautions should always be followed as below:

1. DO NOT attempt to plug in and unplug the unit while feet are in water.
2. Before plugging in or unplugging the unit, be sure the rotary switch is in OFF or "0" position.
3. NEVER immerse the unit completely into water or other liquid for cleaning or other purposes.
4. DO NOT stand up or placing other objects in the unit, use it only when you are in sitting position.
5. DO NOT leave the unit in ON position when the unit are not in use or ready for use.
6. DO NOT operate the unit if it has a damaged cord or plug. If it is not working properly, or has been dropped or damaged in any way whatsoever. Take it to the nearest service centre for examination and repair, or for mechanical adjustment.
7. DO NOT pull the unit by the cord, or use as a handle.
8. CLOSE SUPERVISION is necessary when the unit is used by, or near children, or invalid.
9. DO NOT use while sleeping, or if you feel drowsy.
- 10 USE ONLY on a flat level surface in indoors, never use outdoors.
11. Discontinue use if you feel discomfort, or if any pain or irritation results, never use for swollen or inflamed feet, or if feet have any skin eruptions. First, consult your physician.
12. The appliance has a hot surface. Persons insensitive to heat must be careful when using the appliance.
13. If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard.
14. The appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
15. This appliance may be plugged into AC electrical outlet. (Ordinary household current) Do not use any other outlets.

CARE AND MAINTENANCE

This appliance is for household use purpose, there is no user serviceable parts inside. Do not try to open or change the parts by yourself. Please refer to qualified electrical service or return to the nearest services centre if servicing is needed. This is a health care unit, not a washing basin. Clean feet before use and use clean and soft water. **Deposit will clog air jets and undermine the bubble effect.**

TO CLEAN

After using, simply rinse out with water. Never immerse the entire unit into water. If you wish, you can wipe off all surfaces with a cloth dampened with water and a mild detergent, then rinse and dry.

TO STORE

Be sure the unit are unplugged, cooled and emptied of water. Then store it in a clean and dry place. Before storage just to wore the whole power cord round the wiring holder. Never hang the unit by the power cord.

HOW TO USE THE MASSAGER

The appliance must always be operated with water. Check that the rotary switch is in OFF or "0" position before you connect the appliance to the power.

Do not connect or disconnect the appliance to/from the mains when your feet are in the water or when your hands are wet.

Place the appliance on a flat and suitable surface. The non-skid rubber feet prevent it from slipping. Fill the appliance with water (of desired temperature). Do not fill water over the "MAX" water-level line marked on the unit. This unit is designed to maintain the water temperature at a comfortable level, but not for heating.

Make yourself seat comfortably, place feet on the foot rests and select your desired setting: vibration, bubble, massage with or without warm water.

USING THE CONTROL SWITCH YOU CAN SELECT THE DESIRED APPLICATION

OFF (0)

Turn the unit off.

WARMTH + BUBBLE + MASSAGE + Red light [W + B + M + IR] (Fig. 5) (1)

Turn the rotary switch to position „1“. This function will give you a unique vibration massage and a relaxing footbath in warm water with stimulatory bubble bath. Red light is also activated.

WARMTH + BLOWING [W + A] (2)

Set switch to this mark if you want to massage your feet in warm water at a comfortable temperature.

The unique aroma diffuser provides a personal and convenient aromatherapy treatment. The aroma diffuser provides real aromatherapy, since the aroma is optimal diffused directly by the air blow. You can use a cotton pad and just put a few drops of aromatic oil on the pad and put the pad in the aroma diffuser. Also you can use it as a nail dryer purpose.

MASSAGE + INFRARED [M + IR] (3)

Set switch to this mark if you want to massage your feet in cold water. A foot vibration massage with infrared is also possible with this setting.

When you are through with the massager, just rotate the switch to OFF or "0" position. Then remove your feet and unplug the unit.

REMEMBER NOT TO STAND ON THE UNIT.

Fill with warm water (in desired temperature) into the unit and do not fill water over the "MAX" water-level line on the keel. This unit is designed to maintain the water temperature at a comfortable level, not to raise it.

TO EMPTY, tilt the unit over a sink or tub to enable water to flow out through drain spouts. Rinse out and dry the unit. Never immerse entire unit in water or any other liquid.

Place the appliance on a level surface and the non-skid rubber feet will prevent it from slipping. Be sure the rotary switch is in OFF or "0" position. Then plug the unit into AC outlet (regular household current). Do not stand on the appliance whenever the unit is operated or not. Do not plug the appliance in while your feet are in the water.

SPECIFICATION:

Power Consumption: 230 V/~50 Hz

Output: 90 W

The noise emission of the appliance is less than 70 dB (A)

ECOLOGY – ENVIRONMENTAL PROTECTION

The "crossed-out trash bin" symbol on electrical equipment or packaging indicates that the device cannot be treated as general household waste and should not be disposed of into containers for such waste.

Obsolete or broken-down electrical device should be delivered to special designated collection points, organized by local public administration, whose purpose is to collect recyclable electrical equipment.

This way, each household helps reduce possible negative effects influencing natural environment and allows the reclaiming of materials which the product is made of.

Weight: 2.8kg e

If water leaks from the appliance, the appliance should no longer be used.

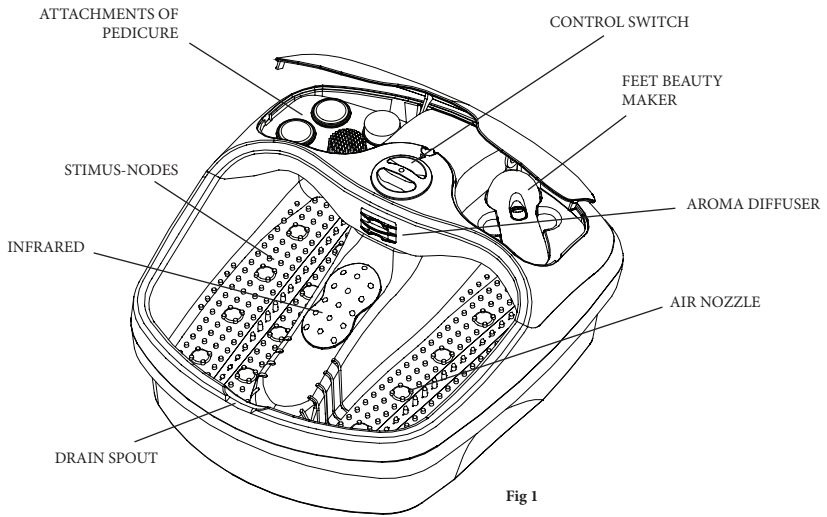
Fig 1

PUT FEW DROPS OF AROMATIC OIL ON THE PAD

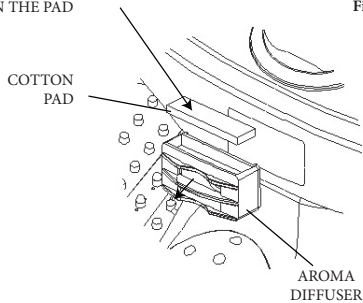
Fig 2

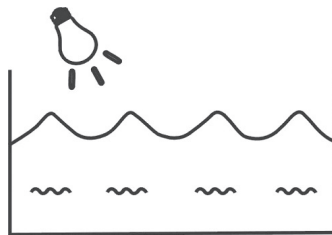
Fig 3

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